


ERGO
TOUCH
GLIDE

IOL INJECTOR KIT

Bulk Purchase

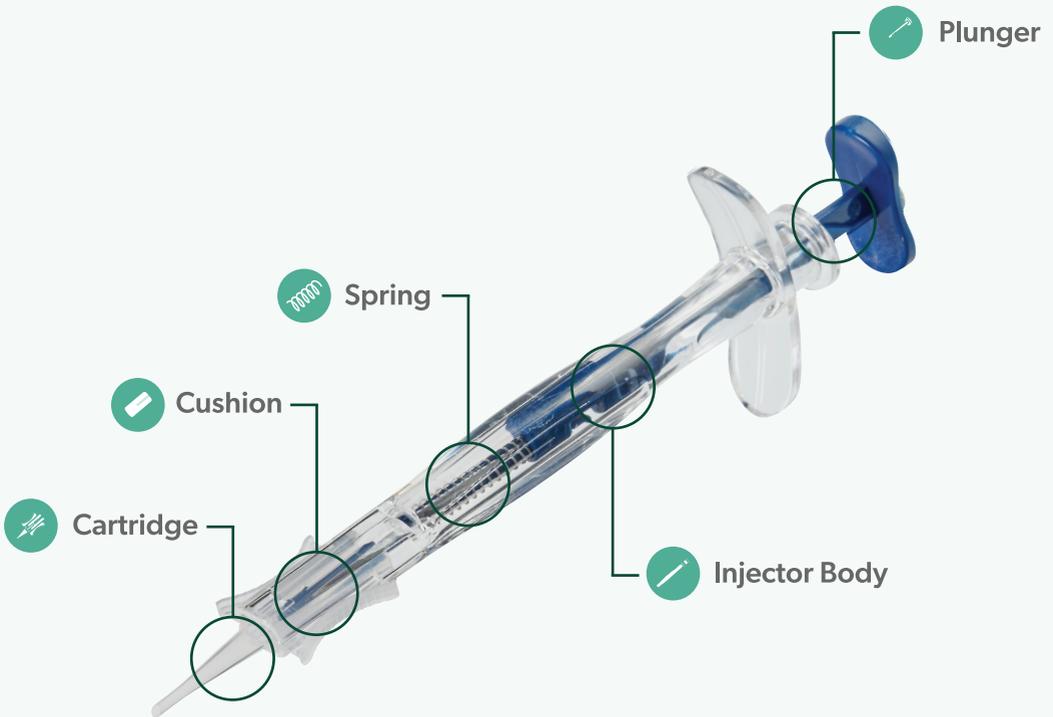

COMPONENTS & DESCRIPTION



Spring

Allows push back for safety & controlled delivery



Plunger

Engineered design & robust material for secure alignment



Injector Body

Designed for both mono & bi-manual injection procedure & allows visibility during delivery process



Cushion

Safety stopper for a controlled start



Cartridge

Cartridge size from 2.2mm to 2.4mm

ASSEMBLY INSTRUCTIONS



Put the **Spring** on the **Plunger**



Insert the **Plunger** into the **Injector Body**



Insert the **Cushion** on top of the **Plunger**



Slide the **Intraocular Lens** into the **Cartridge**



Place the **Cartridge** into the **Injector** and **click** it in place



See step-by-step assembly images

HOW TO LOAD

See the process of loading an IOL into the ErgoTouch GLIDE injector for surgery



Prices available on request and subject to volume

